

LifePort DataStation Case Report

Device Information

Device Information: Kidney Transporter

Serial Number: 1423112

Unit ID: BERN1

Firmware Version: 515

File ID: 39

Data State: Complete

Organ Information

Organ ID: ST-0001-Example

Kidney Side: Right

Blood Type: A

Unofficial Start Time: [REDACTED]

Unofficial Run Time: 04:14:40

Cross Clamp Time & Date:

Cross Clamp Timezone:

Perfusate Used:

Perfusate Lot:

Perfusate Expiration Date:

Cassette Lot #:

Cassette Expiration Date:





Cannula:

Cannula Expiration Date:

Notes:

Organ ID: ST-0001-Example (Right)

This report supplements a handwritten clinical chart. Always refer to official records.

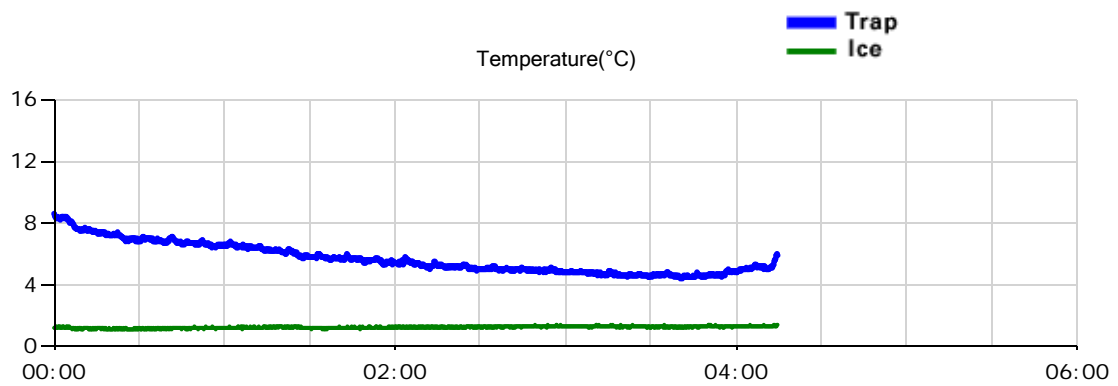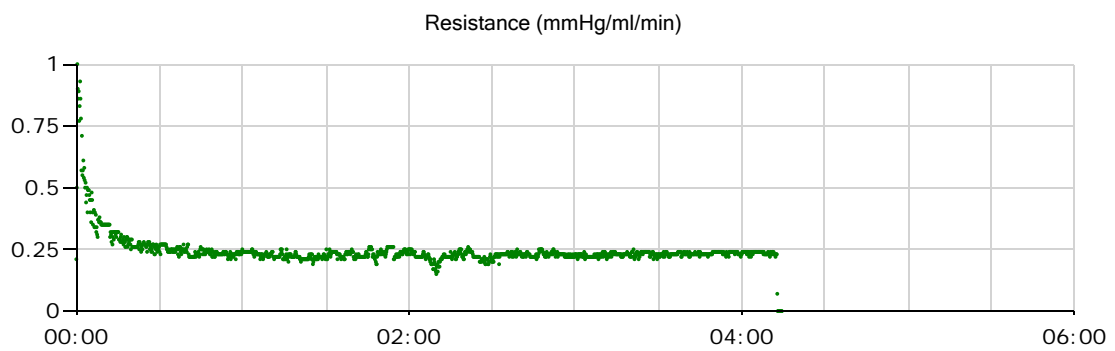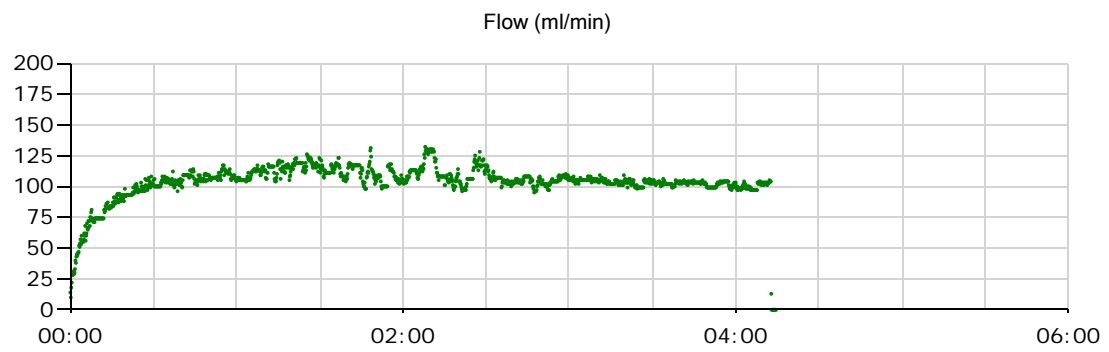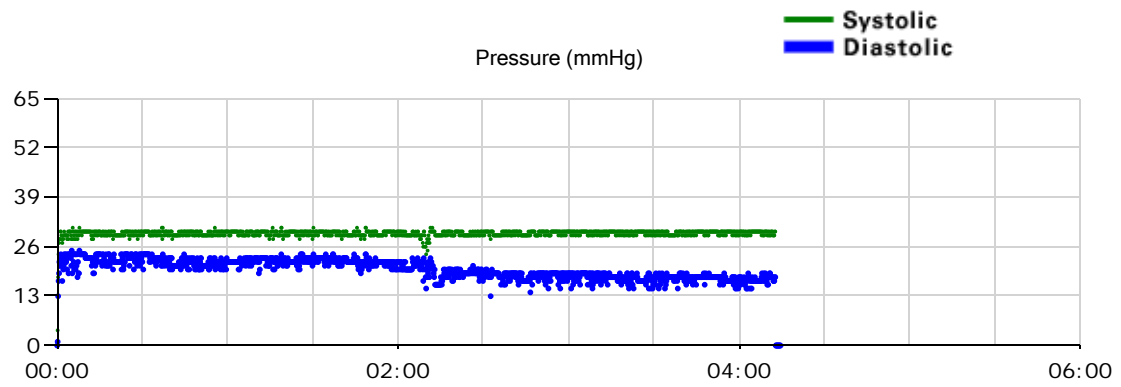

Organ ID: ST-0001-Example (Right)

This report supplements a handwritten clinical chart. Always refer to official records.

Raw Data

Clock Time	Infuse Time	Systolic Press.	Diastolic Press.	Flow	Resist	Ice Temp.	Infuse Temp.
08:44:37	00:00	0	0	14	0.21	1.2	8.6
08:45:37	00:01	28	21	31	0.77	1.2	8.3
08:46:37	00:02	29	23	38	0.71	1.2	8.3
08:59:37	00:15	30	24	84	0.32	1.2	7.4
09:14:37	00:30	29	22	105	0.24	1.2	6.9
09:29:37	00:45	30	23	103	0.25	1.2	6.7
09:44:37	01:00	30	22	105	0.24	1.2	6.6
09:59:37	01:15	29	23	120	0.21	1.2	6.3
10:14:37	01:30	30	23	122	0.22	1.2	5.8
10:29:37	01:45	29	22	114	0.22	1.2	5.7
10:44:37	02:00	30	22	103	0.25	1.2	5.4
10:59:37	02:15	29	18	101	0.24	1.2	5.3
11:14:37	02:30	30	19	115	0.20	1.2	5.0
11:29:37	02:45	29	17	108	0.22	1.3	5.0
11:44:37	03:00	30	18	106	0.23	1.3	4.8
11:59:37	03:15	30	18	102	0.23	1.3	4.7
12:14:37	03:30	30	18	103	0.23	1.2	4.6
12:29:37	03:45	30	18	104	0.24	1.3	4.5
12:44:37	04:00	29	16	103	0.23	1.3	4.8

Event Data

Clock Time	Start Time	Stop Time	Event
12:59:07	04:14:30		Stopped-Pressure Sensor Failure

Organ ID: ST-0001-Example (Right)

This report supplements a handwritten clinical chart. Always refer to official records.

State Data

Clock Time	Start Time	Event
08:44:37	00:00:00	16 (INFUSE)
08:50:17	00:05:40	16 (INFUSE)
08:50:17	00:05:40	128 (LID SENSOR)
08:54:47	00:10:10	48 (INFUSE - Purge Bubbles)
08:54:47	00:10:10	128 (LID SENSOR)
08:54:57	00:10:20	16 (INFUSE)
08:54:57	00:10:20	128 (LID SENSOR)
09:04:57	00:20:20	48 (INFUSE - Purge Bubbles)
09:04:57	00:20:20	128 (LID SENSOR)
09:05:17	00:20:40	16 (INFUSE)
09:05:17	00:20:40	128 (LID SENSOR)
09:15:17	00:30:40	48 (INFUSE - Purge Bubbles)
09:15:17	00:30:40	128 (LID SENSOR)
09:15:27	00:30:50	16 (INFUSE)
09:15:27	00:30:50	128 (LID SENSOR)
09:25:27	00:40:50	48 (INFUSE - Purge Bubbles)
09:25:27	00:40:50	128 (LID SENSOR)
09:25:47	00:41:10	16 (INFUSE)
09:25:47	00:41:10	128 (LID SENSOR)
09:35:47	00:51:10	48 (INFUSE - Purge Bubbles)
09:35:47	00:51:10	128 (LID SENSOR)
09:35:57	00:51:20	16 (INFUSE)
09:35:57	00:51:20	128 (LID SENSOR)
09:45:57	01:01:20	48 (INFUSE - Purge Bubbles)
09:45:57	01:01:20	128 (LID SENSOR)
09:46:17	01:01:40	16 (INFUSE)
09:46:17	01:01:40	128 (LID SENSOR)
09:56:17	01:11:40	48 (INFUSE - Purge Bubbles)
09:56:17	01:11:40	128 (LID SENSOR)
09:56:27	01:11:50	16 (INFUSE)
09:56:27	01:11:50	128 (LID SENSOR)
10:06:27	01:21:50	48 (INFUSE - Purge Bubbles)
10:06:27	01:21:50	128 (LID SENSOR)
10:06:47	01:22:10	16 (INFUSE)
10:06:47	01:22:10	128 (LID SENSOR)

Organ ID: ST-0001-Example (Right)

This report supplements a handwritten clinical chart. Always refer to official records.

10:06:47	01:22:10	128 (LID SENSOR)
10:16:47	01:32:10	48 (INFUSE - Purge Bubbles)
10:16:47	01:32:10	128 (LID SENSOR)
10:16:57	01:32:20	16 (INFUSE)
10:16:57	01:32:20	128 (LID SENSOR)
10:26:57	01:42:20	48 (INFUSE - Purge Bubbles)
10:26:57	01:42:20	128 (LID SENSOR)
10:27:17	01:42:40	16 (INFUSE)
10:27:17	01:42:40	128 (LID SENSOR)
10:37:17	01:52:40	48 (INFUSE - Purge Bubbles)
10:37:17	01:52:40	128 (LID SENSOR)
10:37:27	01:52:50	16 (INFUSE)
10:37:27	01:52:50	128 (LID SENSOR)
10:47:27	02:02:50	48 (INFUSE - Purge Bubbles)
10:47:27	02:02:50	128 (LID SENSOR)
10:47:47	02:03:10	16 (INFUSE)
10:47:47	02:03:10	128 (LID SENSOR)
10:57:47	02:13:10	48 (INFUSE - Purge Bubbles)
10:57:47	02:13:10	128 (LID SENSOR)
10:57:57	02:13:20	16 (INFUSE)
10:57:57	02:13:20	128 (LID SENSOR)
11:07:57	02:23:20	48 (INFUSE - Purge Bubbles)
11:07:57	02:23:20	128 (LID SENSOR)
11:08:17	02:23:40	16 (INFUSE)
11:08:17	02:23:40	128 (LID SENSOR)
11:18:17	02:33:40	48 (INFUSE - Purge Bubbles)
11:18:17	02:33:40	128 (LID SENSOR)
11:18:27	02:33:50	16 (INFUSE)
11:18:27	02:33:50	128 (LID SENSOR)
11:28:27	02:43:50	48 (INFUSE - Purge Bubbles)
11:28:27	02:43:50	128 (LID SENSOR)
11:28:47	02:44:10	16 (INFUSE)
11:28:47	02:44:10	128 (LID SENSOR)
11:38:47	02:54:10	48 (INFUSE - Purge Bubbles)
11:38:47	02:54:10	128 (LID SENSOR)
11:38:57	02:54:20	16 (INFUSE)
11:38:57	02:54:20	128 (LID SENSOR)
11:48:57	03:04:20	48 (INFUSE - Purge Bubbles)
11:48:57	03:04:20	128 (LID SENSOR)

Organ ID: ST-0001-Example (Right)

This report supplements a handwritten clinical chart. Always refer to official records.

11:48:57	03:04:20	128 (LID SENSOR)
11:49:17	03:04:40	16 (INFUSE)
11:49:17	03:04:40	128 (LID SENSOR)
11:59:17	03:14:40	48 (INFUSE - Purge Bubbles)
11:59:17	03:14:40	128 (LID SENSOR)
11:59:27	03:14:50	16 (INFUSE)
11:59:27	03:14:50	128 (LID SENSOR)
12:09:27	03:24:50	48 (INFUSE - Purge Bubbles)
12:09:27	03:24:50	128 (LID SENSOR)
12:09:47	03:25:10	16 (INFUSE)
12:09:47	03:25:10	128 (LID SENSOR)
12:19:47	03:35:10	48 (INFUSE - Purge Bubbles)
12:19:47	03:35:10	128 (LID SENSOR)
12:19:57	03:35:20	16 (INFUSE)
12:19:57	03:35:20	128 (LID SENSOR)
12:29:57	03:45:20	48 (INFUSE - Purge Bubbles)
12:29:57	03:45:20	128 (LID SENSOR)
12:30:17	03:45:40	16 (INFUSE)
12:30:17	03:45:40	128 (LID SENSOR)
12:39:37	03:55:00	16 (INFUSE)
12:40:17	03:55:40	48 (INFUSE - Purge Bubbles)
12:40:27	03:55:50	16 (INFUSE)
12:50:27	04:05:50	48 (INFUSE - Purge Bubbles)
12:50:47	04:06:10	16 (INFUSE)
12:57:27	04:12:50	2 (STOPPED)